



Comfort series

B4w

MANUAL BED



Registration code: SK-CD1



SAIKANG/Medical

Healthier
your life

B4w Manual bed

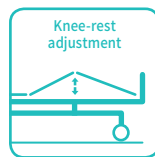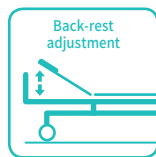

Technical parameters

▪ Size (LxW):	2180x960 mm
▪ Bed height adjustment:	460-710 mm
▪ Mattress platform:	1925x830 mm
▪ Safe working load:	250 kg
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Reverse trendelenburg:	0-12°
▪ Trendelenburg:	0-12°

Technical configuration

☒ Manual crank	4 pcs.
☒ 5" hardshell castors with break	4 pcs.
☒ Foldable side rail with ProErgo-Handle™	1 set
☒ ABS bed ends with safe lock	1 set
☒ IV pole	1 set
☒ IV pole prevision	4 pcs.
☒ Drainage hook	2 pcs.

Optional

☐ Directional wheels function	
☐ Monkey pole	
☐ Aluminium balkan frame	

Design

Patented L-shaped foldable side guardrail: bended tube, solid installing, no shaking, anti-pinch, protect patient's safety.



Hand control

Retractable hand crank, 100,000 times remaining intact. Extreme protection device, sounds when tipped to the max position.



Bed surface

Stamped forming platform, smooth surface, no welding scar, safe working load 250 kg. Optional: ABS platform and resin X-ray translucent platform.



Painting

Baosteel from Fortune Global 500. 11 process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, painting can resist 50 kg impact.



TPR tire castors

TPR tire and save anti-winding hard shell, united forming without bolts. Pass Dynamic test: bearing 120 kg run 30 km and pass obstacles 500 times.



Virgin ABS plastic

Environmental plastic can be sterilized in 100°C. Tensile upto 30MPa.



SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
+86 18021231901 · 24/7
sales@chinasaikang.com
www.saikangmedical.com



Company video



General catalogue