



SK-EX102

Mobile Electric
Suction
Device



SAIKANG / Medical

LEADING BRAND OF CHINA MEDICAL FURNITURE
Stock Code: 870098

TECHNICAL PARAMETERS

- / Power voltage: AC220V \pm 10% 50Hz
- / Negative pressure: ≥ 0.09 MPa(680mmHg)
- / Suction rate: ≥ 20 L/min(35L pump suitable)
- / Range of negative pressure: 0.013MPa \sim 0.09Mpa(680mmHg)
- / Suction bottle: 2500ml x 2
- / Noise: ≤ 55 dB
- / Input power: 400VA
- / Pump structure: Diaphragm type
- / Working type: intermittent loading, continuous operation
- / All - plastic ABS enclosure that is formed in one piece is adopted
- / This device is equipped with a imported membrane pump available with long service life
- / This device has low noise and is maintenance - free
- / This device is equipped with an overflow protection unit to prevent liquid from getting into pump body
- / It's equipped with the hand-operated switch and footstep - operated switch



SWITCH


have manual switch and foot switch, easy to operate



QUALITY CONTROL

OUTWARD

One-time molding of outer ABS case, seamless and durable

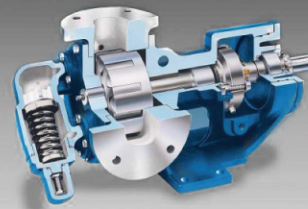

PUMP

Imported membrane pump for long service life

DIAL

Large dial data display for easy reading

Intermittent operation

PRESSURE



ELECT



BOTTLE

5L bottle, suitable for big flow usage requirements

FILTER

Filters improve hygiene standards



PROTECTIVE DEVICE

Overflow protection device prevents liquid from entering the pump



INNOVATION OF SAIKANG

LEADING BRAND OF CHINA MEDICAL FURNITURE



Design Innovation



Research Innovation



Manufacturing Innovation



Maintenance Innovation



Service Innovation

JIANGSU SAIKANG MEDICAL EQUIPMENT CO.,LTD.

24hrs: 0086-18021231901 
www.chinasaikang.com

SAIKANG (China Headquarter)

No.35 Lehong Road,Zhangjiagang City,
Jiangsu Province,China.
TEL: +86-512-58525850
FAX: +86-512-58520713
E-mail: inquiry@jssaikang.com

SAIKANG (USA Branch)

Saikang Medical Technology, Inc
17800 Castleton St Ste 180,
City of Industry, Ca 91748
Tel: +1 949648-2776

SAIKANG (UK Branch)

UK Healthward International Limited
Owner(s) Address-DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +44 2032399738

SAIKANG (East Africa Branch)

Saikang Medical Solutions (East Africa) Ltd.
Viraj Complex Unit 3 Mombasa Road,
P.O Box 14262-00800, Nairobi.
TEL: +254 733-704 690

SAIKANG (Mid East Branch)

Saikang Medical Supplies (Mid East) LLC
X22-S03 England Cluster,
International City, Dubai UAE
TEL: +20 128 8962497



Catalogue



Website

As the circumstance of printing, the photos might have minor difference on color, tactile from real products.
As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.