



A6k

THREE FUNCTIONS ELECTRIC BED

Registration code: SKD-A



SAIKANG Medical

Healthier
your life

A6k

FEATURES AND HIGHLIGHTS




Technical parameters

▪ External size (LxWxH)	2110x960x400-710 mm
▪ Mattress platform	1925x820 mm
▪ Safe working load	200 kg
▪ Back-rest adjustment	0-75° (±10°)
▪ Knee-rest adjustment	0-35° (±10°)
▪ Voltage	100-240 V
▪ Frequency	50/60 Hz
▪ Plug as required (hospital grade)	



Technical configuration

☒ Linear Motor	3pcs
☒ 5" Covered Castors	4pcs
☒ Foldable Side Rail	1set
☒ ABS Bed Ends with Safe Lock	1set
☒ IV Pole	1set
☒ IV Pole Prevision	4pcs
☒ Drainage Hook	4pcs
☒ Mattress Holder	2pcs

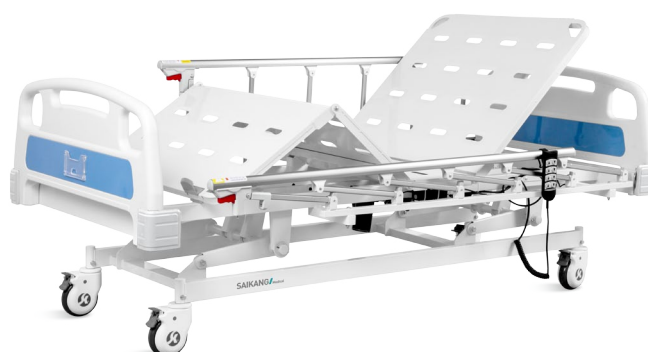

Optional

- | | |
|--|--|
| ☐ Directional Wheel | |
|--|--|



Steels | Painting

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50kg impact.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.

**Bumper Cover**

Reduce the impact of the collision, safe and slinky.

**Electric System**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Design**

Modern design. Use new and strong inserted ABS fire-proof panel avoid off. Easy to be cleaned and installed. Net weight over 3.6kg.

**Platform**

Stamping forming platform, smooth surface no welding scar, bearing 200KG load. Option resin X-ray translucent platform.

**Siderail**

Aluminum alloy side rail, 6-staff, collapsed type with safe lock.

**Castors**

5" covered casters. TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts. Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 time.

**Drainage Hooks**

Laser cutting drainage hook, both sides.

**IV Pole**

Full SS IV pole, bearing 15kg.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service

WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
 +86 18021231901 · 24/7
 sales@chinasaikang.com
 www.saikangmedical.com



Company video



General catalogue