



B6e

THREE FUNCTIONS ELECTRIC BED

Registration code: SKD-C






SAIKANG/Medical

Healthier
your life

B6e

FEATURES AND HIGHLIGHTS

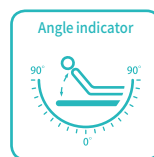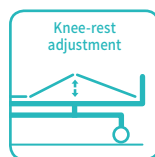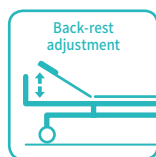

Technical parameters

▪ External size (LxW):	2170x1060 mm
▪ Mattress platform:	1925x900 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	435-745 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Angle indicator:	0-90°
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	



Technical configuration

<input checked="" type="checkbox"/> Linear Motor	3 pcs
<input checked="" type="checkbox"/> Hand Controller	1 pc
<input checked="" type="checkbox"/> 5" Hardshell Castors	4 pcs
<input checked="" type="checkbox"/> ABS Side Rail Extra Height	1 set
<input checked="" type="checkbox"/> ABS Bed Ends with Safe Lock	1 set
<input checked="" type="checkbox"/> IV Pole	1 set
<input checked="" type="checkbox"/> IV Pole Prevision	4 pcs
<input checked="" type="checkbox"/> Drainage Hook	2 pcs




Optional

<input type="checkbox"/> Manual CPR	
<input type="checkbox"/> Electric CPR	
<input type="checkbox"/> Back-up Battery	
<input type="checkbox"/> X-ray Bed Surface	
<input type="checkbox"/> Monkey Pole	
<input type="checkbox"/> Balkan Traction Frame	
<input type="checkbox"/> Directional Wheel	



Design

E-shaped large guardrail, with angled display and built-in controller, meets safety standards, providing comforting feeling for patients.

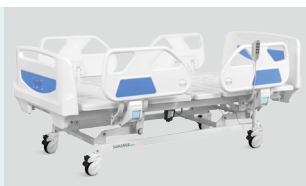**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding. Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**


Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Virgin ABS**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**Angle Indicator**

Patented angle indicator.

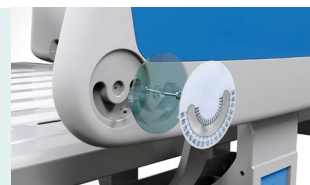**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

**Side Rail**

Full ABS material, tuck away side rail with built-in controller. From top of the guardrail to bed surface is about 345mm, more height provides better protection for patient. Bottom can endure 100,000 times press, buffer, built-in angle indicator.

**Head & Foot Board**

Full ABS material, Bottom can endure 100,000 times press. Head and foot board can be reversible.

**Wheels**

TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts. Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 times.

**IV Pole**

Full stainless steel IV pole, bearing 15 kg.

**Drainage Hooks**

Mobile plastic drainage hooks.





CONSTANT
INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE
BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ sales@chinasaikang.com
🌐 www.saikangmedical.com



Company video



General catalogue