



B8e

MULTIFUNCTIONAL ICU ELECTRIC BED

Registration code: SKD-C

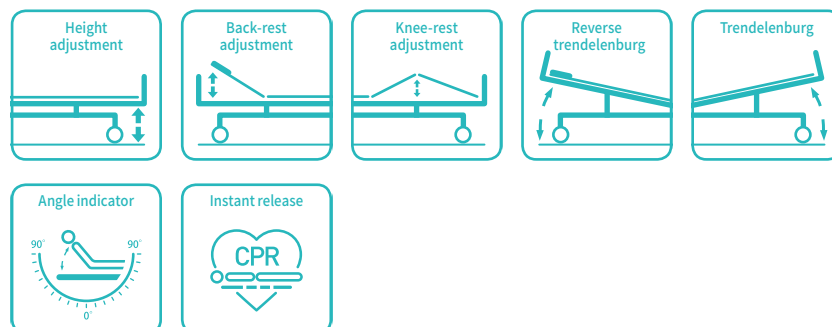


SAIKANG/Medical

Healthier
your life

B8e

FEATURES AND HIGHLIGHTS



Technical parameters

▪ External size (LxW):	2170x1065 mm
▪ Mattress platform:	1925x900 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	435-725 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Reverse Trendelenburg:	0-12°
▪ Trendelenburg:	0-12°
▪ Angle indicator:	0-90°
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	



Technical configuration

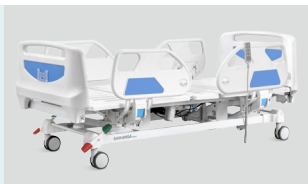
<input checked="" type="checkbox"/> Linear Motor	4 pcs
<input checked="" type="checkbox"/> 5" Double Side Castors, Central Locking	4 pcs
<input checked="" type="checkbox"/> ABS Side Rail Extra Height	1 set
<input checked="" type="checkbox"/> ABS Bed Ends with Integrated Control Panel	1 set
<input checked="" type="checkbox"/> IV Pole	1 set
<input checked="" type="checkbox"/> IV Pole Prevision	4 pcs
<input checked="" type="checkbox"/> Drainage Hook	2 pcs
<input checked="" type="checkbox"/> Central Locking Pedal	2 pcs
<input checked="" type="checkbox"/> Electric CPR	

Optional

<input type="checkbox"/> Manual CPR
<input type="checkbox"/> Weighing System
<input type="checkbox"/> Back-up Battery
<input type="checkbox"/> X-ray Bed Surface
<input type="checkbox"/> Monkey Pole
<input type="checkbox"/> Balkan Traction Frame
<input type="checkbox"/> Directional Wheel

Design

E-shaped large guardrail, with angled display and built-in controller, meets safety standards, providing comforting feeling for patients.

**Steels | Painting**

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**

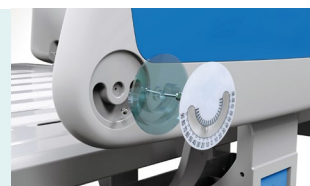
Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Virgin ABS**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**Angle Indicator**

Patented angle indicator.

**Wheels**

Secure® Castors
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**Durable Brakes**

Break spring keep qualified after 10,000 times uninteruped break.

**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

**Side Rail**

Full ABS material, tuck away side rail with built-in controller. From top of the guardrail to bed surface is about 345mm, more height provides better protection for patient. Bottom can endure 100,000times press, buffer, built-in angle indicator.

**Head & Foot Board**

Full ABS material, Bottom can endure 100,000times press
Head and foot board can be reversible.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV Pole**

Full stainless steel IV pole, bearing 15 kg.



SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

+86 18021231901 · 24/7

sales@chinasaikang.com

www.saikangmedical.com



Company video



General catalogue