



Grace Care series

# D2w

## TWO FUNCTIONS MANUAL BED

Registration code: SK-CD1




**SAIKANG** / Medical

Healthier  
your life

## D2w

### FEATURES AND HIGHLIGHTS

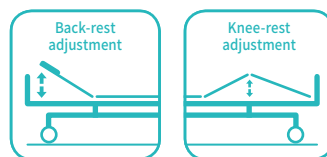

#### Technical parameters

|                         |                  |
|-------------------------|------------------|
| ▪ External size (LxWxH) | 2245x1000x510 mm |
| ▪ Mattress platform     | 1925x830 mm      |
| ▪ Safe working load     | 250 kg           |
| ▪ Back-rest adjustment  | 0-75° (±10°)     |
| ▪ Knee-rest adjustment  | 0-35° (±10°)     |



#### Technical configuration

|   |      |
|---|------|
| <input checked="" type="checkbox"/> Manual Cranks                           | 2pcs |
| <input checked="" type="checkbox"/> 5" Hardshell Castors                    | 4pcs |
| <input checked="" type="checkbox"/> Foldable Side Rail with ProErgo-Handle™ | 1set |
| <input checked="" type="checkbox"/> ABS Grace Bed Ends with Safe Lock       | 1set |
| <input checked="" type="checkbox"/> IV Pole                                 | 1set |
| <input checked="" type="checkbox"/> IV Pole Prevision                       | 4pcs |
| <input checked="" type="checkbox"/> Drainage Hook                           | 2pcs |
| <input checked="" type="checkbox"/> Shoes Holder                            | 1set |



#### Optional

|  |  |
|--|--|
| <input type="checkbox"/> Monkey Pole           |  |
| <input type="checkbox"/> Balkan Traction Frame |  |
| <input type="checkbox"/> Directional Wheel     |  |



**Design**

Extra long and height head foot board, with curve surrounded design, provide comforting feeling for patients.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.

**Steels | Painting**

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding Solid bed base made of 5.0 mm stamped rotating plate.

**Hand Control**

Retractable hand crank, 100,000 times remaining intact  
Extreme protection device, sounds when tipped to the max position,  
Strong pallets to increase the bed height adjustment.

**Platform**

Stamping forming platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform.

**Guardrail**

L-shaped foldable guardrail (Patent) with ProErgo-Handle™, patented parts. with red and green mark to indicated if the siderail is locked or not. Bended tube, solid installing, 0 shaking, anti-pinch, protect patient's safety.

**Head & Foot Board**

Head and foot board can be reversible, extra height design, size: 1000\*563mm; ABS Grace Bed Ends with Safe Lock, automatic lock without any manual operation on lock, on/off removable mark on hook lock.

**Castors**

5" hardshell casters  
TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts.  
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 time.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV Pole**

Full SS IV pole, bearing 15kg.

**Shoes Holder**

Shoes holder: 800\*400mm, adjustable





## **CONSTANT INNOVATION**



Design



Research



Manufacturing



Maintenance



Service

## **WORLDWIDE BRANCHES**

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
 +86 18021231901 · 24/7  
 sales@chinasaikang.com  
 www.saikangmedical.com



Company video



General catalogue