




Grace Care series

# D3d

## THREE FUNCTIONS MANUAL BED

Registration code: SK-CD1

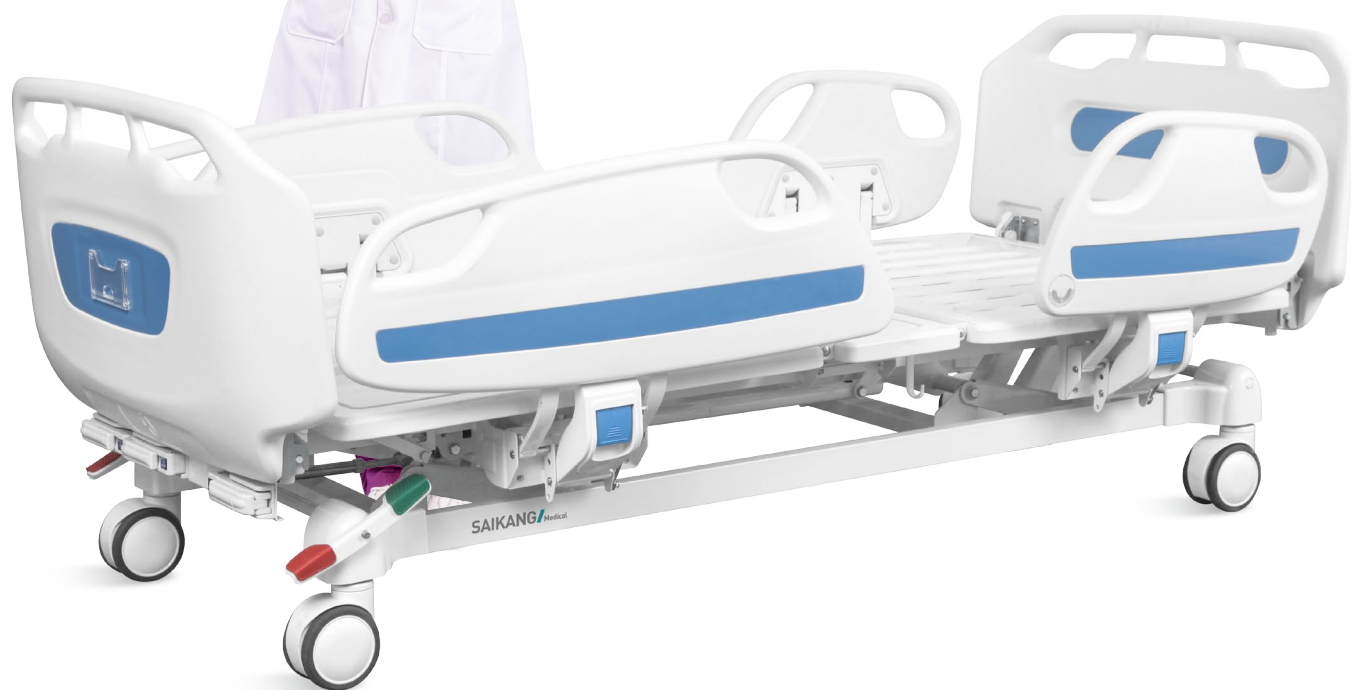

**SAIKANG** Medical

Healthier  
your life


# D3d

## FEATURES AND HIGHLIGHTS



### Technical parameters

|                         |                      |
|-------------------------|----------------------|
| ▪ External size (LxWxH) | 2245x1090x435-745 mm |
| ▪ Mattress platform     | 1925x900 mm          |
| ▪ Safe working load     | 250 kg               |
| ▪ Back-rest adjustment  | 0-75° (±10°)         |
| ▪ Knee-rest adjustment  | 0-35° (±10°)         |

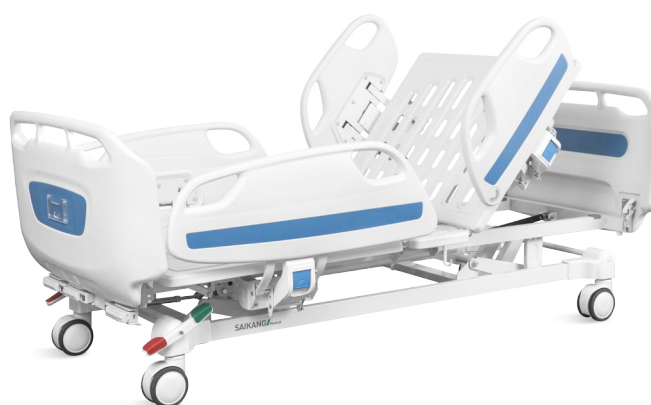

### Technical configuration

|                                                                             |      |
|-----------------------------------------------------------------------------|------|
| <input checked="" type="checkbox"/> Manual Cranks                           | 3pcs |
| <input checked="" type="checkbox"/> 5" Double Side Castors, Central Locking | 4pcs |
| <input checked="" type="checkbox"/> ABS Grace Side Rail Extra Height        | 1set |
| <input checked="" type="checkbox"/> ABS Grace Bed Ends with Safe Lock       | 1set |
| <input checked="" type="checkbox"/> Central Locking Pedal                   | 2pcs |
| <input checked="" type="checkbox"/> IV Pole                                 | 1set |
| <input checked="" type="checkbox"/> IV Pole Prevision                       | 4pcs |
| <input checked="" type="checkbox"/> Drainage Hook                           | 2pcs |




### Optional

|                                                |  |
|------------------------------------------------|--|
| <input type="checkbox"/> Monkey Pole           |  |
| <input type="checkbox"/> Balkan Traction Frame |  |
| <input type="checkbox"/> Directional Wheel     |  |



**Design**

Extra long and height head foot board, with curve surrounded design, provide comforting feeling for patients.

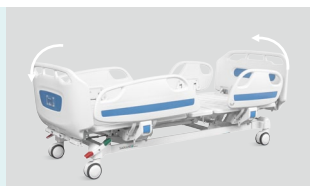**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.

**Steels | Painting**

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding Solid bed base made of 5.0 mm stamped rotating plate.

**Hand Control**

Retractable hand crank, 100,000 times remaining intact  
Extreme protection device, sounds when tipped to the max position,  
Strong pallets to increase the bed height adjustment (435-745mm).

**Platform**

Stamping forming platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform.

**Guardrail**

Tuck away side rail with angle indicator.

**Head & Foot Board**

Head and foot board can be reversible, extra height design, size: 1000\*563mm; ABS Grace Bed Ends with Safe Lock, automatic lock without any manual operation on lock, on/off removable mark on hook lock.

**Wheels**


TPR tire no worn out after running 30 km, save anti-winding hard shell, united forming without bolts. Pass dynamic test: Bearing 120 kg, run 30 km and pass obstacles 500 time.

**Durable Brakes**

Break spring keep qualified after 10,000 times uninterrupted break.

**Angle Indicator**

Patented angle indicator.

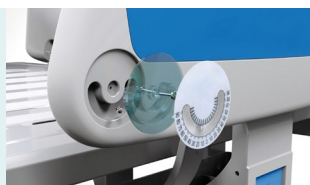**Drainage Hooks**

Mobile plastic drainage hooks.

**IV Pole**

Full SS IV pole, bearing 15kg.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎ +86 18021231901 · 🕒 24/7  
✉ sales@chinasaikang.com  
🌐 www.saikangmedical.com



Company video



General catalogue