



Grace Care series

# D3w

## THREE FUNCTIONS MANUAL BED

Registration code: SK-CD1



**SAIKANG** Medical

Healthier  
your life

## D3w

### FEATURES AND HIGHLIGHTS




#### Technical parameters

▪ External size (LxWxH)	2245x1000x435-745 mm
▪ Mattress platform	1925x830 mm
▪ Safe working load	250 kg
▪ Back-rest adjustment	0-75° ( $\pm 10^\circ$ )
▪ Knee-rest adjustment	0-35° ( $\pm 10^\circ$ )



#### Technical configuration

<input checked="" type="checkbox"/> Manual Cranks	3pcs
<input checked="" type="checkbox"/> 5" Hardshell Castors	4pcs
<input checked="" type="checkbox"/> Foldable Side Rail with ProErgo-Handle™	1set
<input checked="" type="checkbox"/> ABS Grace Bed Ends with Safe Lock	1set
<input checked="" type="checkbox"/> IV Pole	1set
<input checked="" type="checkbox"/> IV Pole Prevision	4pcs
<input checked="" type="checkbox"/> Drainage Hook	2pcs

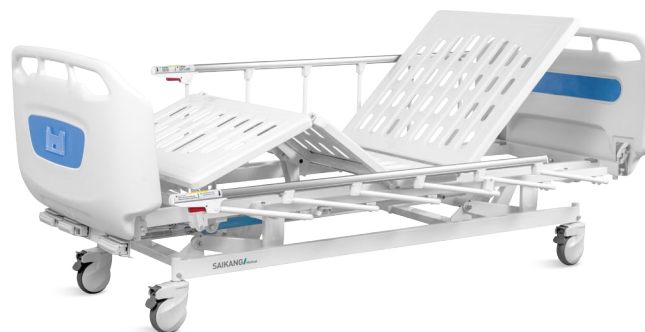

#### Optional

<input type="checkbox"/> Monkey Pole	
<input type="checkbox"/> Balkan Traction Frame	
<input type="checkbox"/> Directional Wheel	



### Design

Extra long and height head foot board, with curve surrounded design, provide comforting feeling for patients.



### Plastic

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.



### Steels | Painting

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.



### Welding | Chassis

Panasonic robotic ensure, 360° full smooth wedding Solid bed base made of 5.0 mm stamped rotating plate.




### Hand Control

Retractable hand crank, 100,000 times remaining intact  
Extreme protection device, sounds when tipped to the max position,  
Strong pallets to increase the bed height adjustment (435-745mm).



### Platform

Stamping forming platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform.

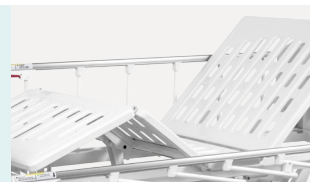

### Guardrail

L-shaped foldable guardrail (Patent) with ProErgo-Handle™, patented parts, with red and green mark to indicated if the siderail is locked or not. Bended tube, solid installing, 0 shaking, anti-pinch, protect patient's safety.



### Head & Foot Board

Head and foot board can be reversable, extra height design, size: 1000\*563mm; ABS Grace Bed Ends with Safe Lock, automatic lock without any manual operation on lock, on/off removable mark on hook lock.



### Castors

5" hardshell casters  
TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts.  
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 time.



### Drainage Hooks

Mobile plastic drainage hooks.



### IV Pole

Full SS IV pole, bearing 15kg.





CONSTANT  
INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE  
BRANCHES

#### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

#### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

#### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

#### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

#### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SAIKANG** Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎ +86 18021231901 · 🕒 24/7  
✉ sales@chinasaikang.com  
🌐 www.saikangmedical.com



Company video



General catalogue