

Grace Care series

D8d


MULTIFUNCTIONAL ELECTRIC BED

Registration code: SKD-C



D8d

FEATURES AND HIGHLIGHTS

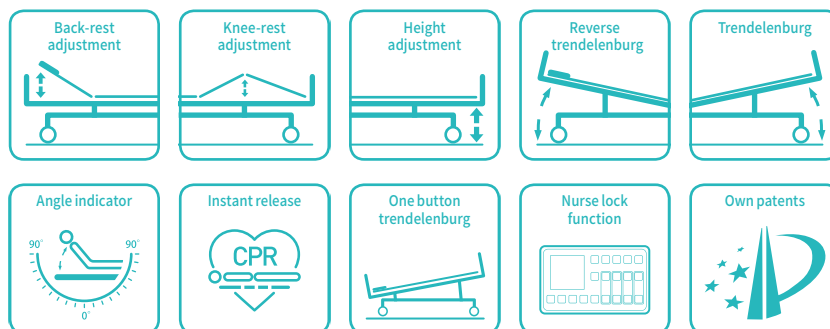

Technical parameters

▪ External size (LxW):	2280x1090 mm
▪ Mattress platform:	1925x900 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	435-725 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Reverse Trendelenburg:	0-12°
▪ Trendelenburg:	0-12°
▪ Angle indicator:	0-90°
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	



Technical configuration

<input checked="" type="checkbox"/> Manual CPR	
<input checked="" type="checkbox"/> One button stop	
<input checked="" type="checkbox"/> One button trendelenburg	
<input checked="" type="checkbox"/> Nurse lock system	
<input checked="" type="checkbox"/> Angle indicator for back-rest and Trendelenburg position	4 pcs.
<input checked="" type="checkbox"/> Linear motor	4 pcs.
<input checked="" type="checkbox"/> 5" double side castors with central locking system	4 pcs.
<input checked="" type="checkbox"/> ABS side rail with integrated control panel	1 set
<input checked="" type="checkbox"/> ABS bed ends with integrated LCD control panel	1 set
<input checked="" type="checkbox"/> φ85 mm wall bumper wheel	4 pcs.
<input checked="" type="checkbox"/> Trapeze socket	2 pcs.
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> IV pole prevision	4 pcs.
<input checked="" type="checkbox"/> Drainage hook	2 pcs.
<input checked="" type="checkbox"/> Central locking pedal	2 pcs.

Optional

<input type="checkbox"/> Electrical CPR	
<input type="checkbox"/> Weighing system with LCD Screen	
<input type="checkbox"/> Bed extension	
<input type="checkbox"/> Back-up battery	
<input type="checkbox"/> Cardiac position	
<input type="checkbox"/> ABS bed surface	
<input type="checkbox"/> X-ray bed surface	
<input type="checkbox"/> Monkey pole	
<input type="checkbox"/> Balkan traction frame	
<input type="checkbox"/> Directional wheel	
<input type="checkbox"/> Basic foam mattress (12 cm thickness)	
<input type="checkbox"/> Anti-decubitus foam mattress (14 cm thickness)	
<input type="checkbox"/> Vertical oxygen holder	
<input type="checkbox"/> Single arm lifting pole	
<input type="checkbox"/> Utility shelf	

Head and Foot Board

Extra long and height head foot board, with curve surrounded design, provide comforting feeling for patients.

**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Steels | Painting**

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Virgin ABS**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**LCD**

Unique design head foot board with built-in controller. 4.3 inch LCD display, shows moving angle. Button can endure 100,000 times press.

**Side Rail**

Tuck away side rail with buffer, built-in controller and angle indicator.

**Wheels**


TPR tire no worn out after running 30 km, save anti-winding hard shell, united forming without bolts. Pass dynamic test: Bearing 120 kg, run 30 km and pass obstacles 500 time.

**Durable Brakes**

Break spring keep qualified after 10,000 times uninterrupted break.

**Angle Indicator**

Patented angle indicator.

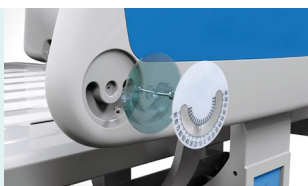**Bumper Wheels**

Patented bumper wheel detachable, max. protect the bed and patient during transportation.

**Safety**

IEC60601 Safety Standard.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV Pole**

Full stainless steel IV pole, bearing 15 kg.

**Optional**

Battery: 18 times movement with 100 kg bearing, 130 times movement without bearing.
Nurse control board.
Storage rack under foot board drawer.
Store control board or linens.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service

WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
 +86 18021231901 · 24/7
 sales@chinasaikang.com
 www.saikangmedical.com



Company video



General catalogue