



X9x

MULTIFUNCTIONAL ELECTRIC BED

Registration code: SKD-C



SAIKANG/Medical

Healthier
your life

X9x


FEATURES AND HIGHLIGHTS

Technical parameters

▪ External size (LxW):	(2040-2240)x1020 mm
▪ Mattress platform:	(1820-2020)x820 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	500-950 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-20° (±10°)
▪ Reverse Trendelenburg:	0-18°
▪ Trendelenburg:	0-18°
▪ Angle indicator:	0-90°
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	

Technical configuration

<input checked="" type="checkbox"/> Linear Motor	6 pcs
<input checked="" type="checkbox"/> 6" Double Side Castors, Central Locking	4 pcs
<input checked="" type="checkbox"/> ABS Side Rail with Integrated Control Panel	1 set
<input checked="" type="checkbox"/> ABS Bed Ends One Step Removable	1 set
<input checked="" type="checkbox"/> Ø75mm Wall Bumper Wheel	4 pcs
<input checked="" type="checkbox"/> IV Pole	1 set
<input checked="" type="checkbox"/> IV Pole Prevision	4 pcs
<input checked="" type="checkbox"/> Drainage Hook	4 pcs
<input checked="" type="checkbox"/> 5th Directional Wheel Pedal	2 pcs
<input checked="" type="checkbox"/> Central Locking Pedal	2 pcs
<input checked="" type="checkbox"/> X-ray Translucent Mattress Platform	1 set
<input checked="" type="checkbox"/> Back-up Battery	1 set
<input checked="" type="checkbox"/> Night Light	1 set
<input checked="" type="checkbox"/> Bed Push Handle for Emergency Transfer	1 set
<input checked="" type="checkbox"/> Manual CPR	
<input checked="" type="checkbox"/> Electric CPR	

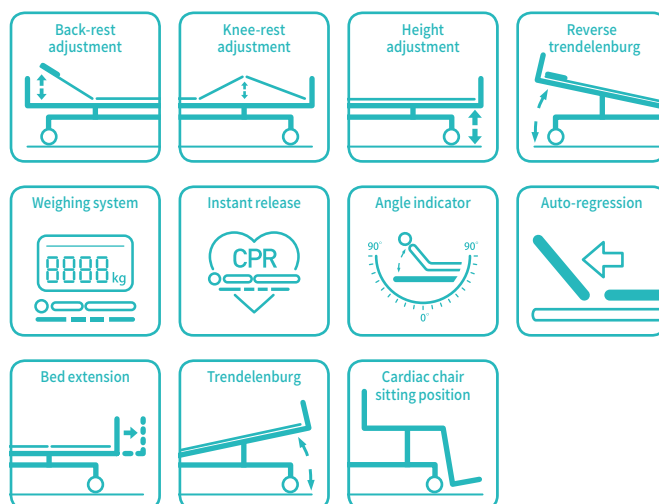

Optional

<input type="checkbox"/> Weighing System With LCD Screen	
<input type="checkbox"/> Data Storage System	
<input type="checkbox"/> Monkey Pole	

Design

Extra Long and height head foot board, with curve surrounded design, provide conforting feeling for patients.

**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

**Welding | Safety**

Panasonic Robotic ensure, 360° full smooth wedding. IEC606011 Safty Standard.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**Electric Motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**7 Inch Touch Screen**

Replacing the original button screen, more convenient operation. There are two language options, chinese and english. Screen can store the patient and doctor's information, users manuals, weight record

**Brake**

Break spring keep qualified after 10,000times uninteruped break.

**Wheels**

TPR tire no worn out after running 30 km, save anti-winding hard shell, united forming without bolts. Pass dynamic test: Bearing 120 kg, run 30 km and pass obstacles 500 time.

**Bed Surface**

X-Ray Translucent Mattress Platform.

**Side Rail**

Tuck away side rail with buffer Built-in controller and angle indicator.

**ABS Bed Ends**

ABS bed ends One Step Removable, meet the first-aid requirements of high-end wards and ICU room.

**Bumper Wheels**

Bumper wheel , $\phi 74$ mm, max. protect the bed and patient during transportation.

**Battery**

Battery - 18times movement with 100kg bearing, 130times movement without bearing.

**IV Pole**

Full stainless steel IV pole, bearing 15 kg.

**Bottom**

The shelf at the bottom of the bed can hold patients items, like shoes and something others.

**Lifting Pole**

Lifting pole - Bearing 765N force.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ sales@chinasaikang.com
🌐 www.saikangmedical.com



Company video



General catalogue